


MAEA Ultrasonic Skin Scrubber



User Manual

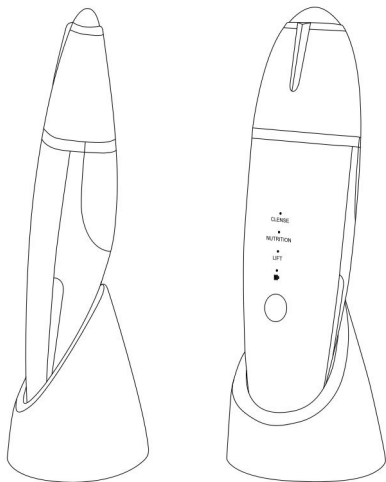

1. People who use implanted medical devices such as pacemakers
2. People who use electrical medical equipment for maintaining life, such as artificial hearts and artificial lungs
3. People using wearable electrical medical devices such as electrocardiogram scanners

Warning: For fear of physical problems, people who meet the following content, please consult a doctor before use



1. People with allergic constitution
2. People with poor blood circulation
3. People who have sensory disturbances such as contagious dermatitis or abnormal perception
4. Those who are taking medication
5. Patients with malignant tumors
6. People with brain abnormalities such as epilepsy
7. People with cognitive impairment due to diabetes, etc.
8. Patients with fever over 38 °C
9. People who need quiet therapy
10. Patients with severe pain such as bone marrow fractures and sprains
11. People with low bone density such as osteoporosis
12. People with abnormal blood pressure
13. People with heart disease
14. Patients under treatment
15. People who feel uncomfortable during use

If you feel unwell even if you do not fit the above conditions, please consult your physician.

Warning: Do not use this product on the following areas. May cause skin or body damage.



1. Eyes, mouth and other mucous membranes
2. Eyeball and upper eyelid
3. Parts undergoing plastic surgery
4. wound
5. Cosmetics and other parts that cause dermatitis
6. Body and other parts other than the face
7. Sites of inflammation due to acne or discharge
8. Metal or plastic embedded in the skin

Do not use in places other than those specified in this instruction manual

Before Nursing

- Confirmation of ultrasonic generation is performed. Moisten the skin with water and slide slowly against the spatula tip. If the water mist flows from the surface of the skin, ultrasonic waves are generated.
- Please wash your face to keep your skin in a beautiful condition.
- Please confirm whether there is any inflammation such as acne and small bumps. Please do not use on areas with inflammation.
- In order to prevent skin problems, an allergy test will be performed in advance. Apply beauty essence or lotion on the inner arm, and check whether there is any abnormality on the skin under each mode. In case of abnormality, please discontinue use and consult a professional doctor.


About Maintenance And Storage

- Please do not store it in a place with high humidity or moisture, or in the accessories of a microwave oven.
After cleaning the main unit, close the cover for storage.
- When dirt or moisture remains on the blade tip, wipe it off with a towel or paper towel.
- When not in use for a long time, in order to prevent aging, please avoid high temperature and humid places, and store the product in a fully charged state.

About Charging

- After installing the main unit, plug the included USB cable into the charging stand and connect to a home computer or AC adapter for charging.
- When charging from an outlet, prepare an AC adapter separately.
- The charging time is about 6 hours. The battery light blinks red during charging. The battery indicator lights up green when charging is complete
- A set of 1 minutes, fully charged can use about 3-4 sets.

Sample Diagram

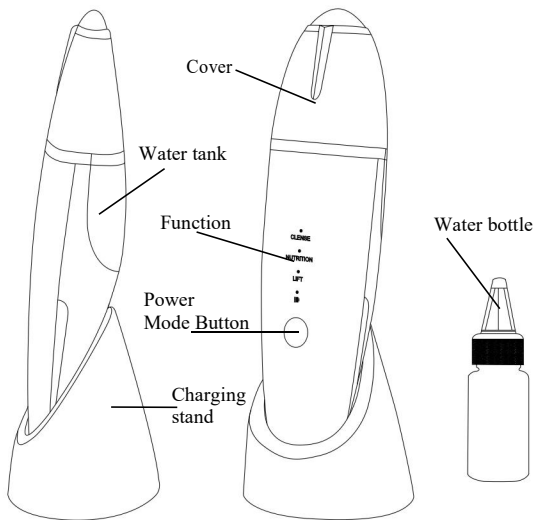

Charging Method

When the power becomes low,
The battery indicator light flashes red,
power off immediately



Connect the USB cable to the charging stand, install the main unit on the charging stand to charge.

The battery indicator light will blink red when charging.

When charging is complete, the indicator light will light up green.



The light will flash red when charging



After charging is complete,
the light turns green

Modes and Button Operations

① Standby mode:

Press the button once to enter standby mode. The "CLEASE", "MUTRITION" and "LIFT" indicators light up in this mode



② Mist mode: Spray

Press the button twice to enter mist mode. In this mode, the "CLEASE" light flashes and water mist is sprayed from the main body.

*If the main body does not tilt downward, the mist will not be sprayed.



③ Water stripping mode: Spray+Vibration

Press the button three times to enter water peel mode.

In this mode, the "CLEASE" indicator lights up. Water mist is sprayed from the main body, and the skin surface is cleaned by ultrasonic vibration.



④ Ion export mode: Iron export+Vibration

Press the button 4 times to enter the ion export mode.

In this mode the "MUTRITION" indicator lights up. The positive micro-current is applied to absorb negatively charged pore dirt



⑤ Ion introduction mode:

Ion introduction mode: Press the button 5 times to enter the ion introduction mode. In this mode, the "NUTRITION" light is on. Stimulate the surface of the skin with micro-current, so that the beauty ingredients can easily penetrate into the stratum corneum of the skin.



⑥ Standby mode:


Press the button 6 times to enter standby mode again.

In this mode, the lights of "CLEASE", "NUTRITION" and "LIFT" are on. After pressing the button, the process will change to the same mode.



⑦ Power off mode:

Power off: long press the button, the power is off. The device automatically powers off after 15 minutes.

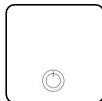

There will be a "di" sound time reminder every 3 minutes, and when the automatic shutdown time is up in 15 minutes, it will automatically shut down after there "didi" alarms.

Blue LED function:

When using each mode, the blue LED light emits blue light from the head. Adjust the texture of the skin and keep it healthy.

Instructions

Fog Mode:


Before use, remove the cover on the outer head of the main unit and use the water bottle to give water. If the subject is tilted downward, it will spray fog for a certain period of time. Do not use substances other than water such as lotion, gel, etc., so as not to cause malfunction.

CLEANS Mode:

Use the mist function to wet the entire face.

As shown in the picture below, face the inner side of the face, keep the front end of the head roughly vertical to the skin, and slide slowly for about 3 minutes.


slide the head towards
the inside of the face

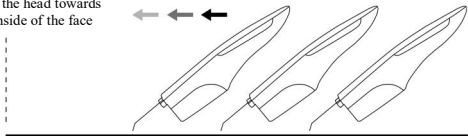

Ion Export Mode:

Before use, apply water, lotion, nutrient solution, etc. to the entire face. As shown in the picture below, move toward the inside of the face, let the front of the head be parallel to the skin, and slide slowly for about 3 minutes.


slide the head towards
the inside of the face

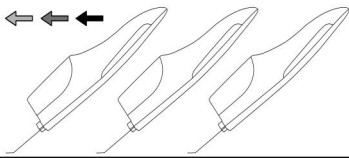

Ion Import Mode:

Before use, apply beauty essence, gel, and lotion to the entire face. As shown in the picture below, toward the outside of the face, let the front of the head be parallel to the skin, and slide slowly for about 3 minutes.

slide the head towards
the outside of the face

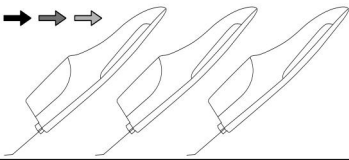

Application area



Facial maintenance

Swipe outward from the center of the face. Move like a glide from the corners of the mouth and the sides of the nose, to the front of the ears, and even from above the eyebrows to the hairline of the hair.



Care of the bridge of the nose and around the nose

Let it slide along the nose from the side to the center of the nose. Then, move up from the front of the nose.




Care around the mouth

Move like a swipe from the jaw to the corners of the mouth, and even from the sides of the mouth to the cheeks.

According to the state of your skin, focus on the parts you care about.

About the water mist function

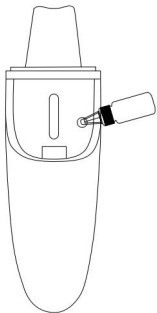

Water mist function

Remove the cover from the outside of the host.

Use the supplied water bottle to supply water from the water supply port.

Replenish each time the water is used up.

Products Include

1. One piece of MAEA Scrubber machine
2. One USB interface Type-C charging cable
3. One plastic water bottle
4. One Manual

Troubleshooting

Content	Possible reason	Solution
Power cannot be turned on	Insufficient charge	please charge
	Battery consumption	Please consult our store
	Failure caused by moisture entering the interior	Please consult our shop
It's too much peeling	Can't change mode	Please refer to P5 mode and button operation
	head slipping fast	please slide slowly
	The correct direction of the slippery head is out of shape	Refer to the mode description, please adjust the direction of the slider
	Not enough moisture on the face	Please apply more water
	battery consumption	Please consult our shop
Can't feel the vibration	Can't change mode	Please refer to P5 mode and button operation
	pushing the slider hard	relax, please slide your head gently
	Insufficient charge	please charge
	battery consumption	Please consult our shop
Dry, burning skin	Proper orientation of the nozzle for dry skin	Please apply more water
	The correct direction of the slippery head is out of shape	Refer to the mode description, please adjust the direction of the slider

If you have any questions about the above matters, please feel free to consult our store.

Free replacement in case of initial defect or failure.

Q: Can it be used on other than the face (body, etc.)?

A: Yes! It's great for use on your hands, feet and body to address the same issues with skin tone or dryness that you would experience on your face.

Q: Is it better to apply water or gel?

A: If you use it on dry skin, it may add burden to the skin. In addition, it is recommended to use it with water or gel in order to transmit ultrasonic waves easily.

Q: How many times can it be used in one day?

A: Within 15 minutes for 1 set, about 2-3 sets can be used as a reference. Please do not use it multiple times in the same place, and do not use it frequently in a short period of time.

Q: Is it waterproof?

A: This product is not waterproof.

Q: Can I return or exchange?

A: In the case of initial defective products, etc., you can return and replace them for you.

Q: How long does this take to improve the skin?

A: If you use it continuously for 2-3 months, you will find that your skin will improve.

Product name:MAEA Skin Scrubber Machine

Host Weight: 130g

Host+charging base:199g

Size: 185*53*36mm

Charging voltage: DC5V, 1A

Power: 2.2W

Battery capacity: 800mAh